


memmert

trust the best

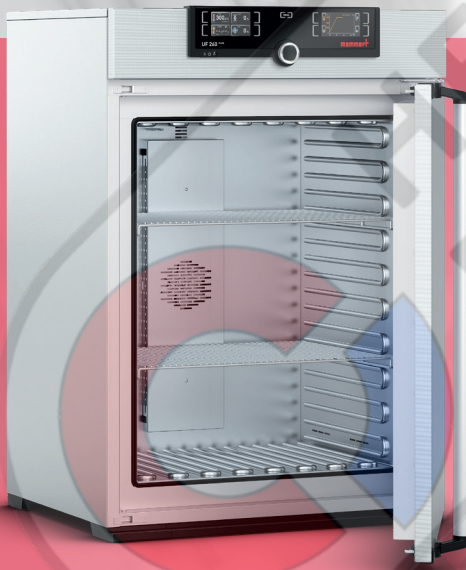

Universal Oven

UN | UF

**Intuitive.
Flexible.
Powerful.**

www.memmert.com

Universal Oven U

IDEAL FOR DRYING, HEATING, TEMPERATURE CONTROL, AGING, BURN-IN, TESTING AND CURING VARIOUS COMPONENTS FROM ALL AREAS OF INDUSTRY AND QUALITY ASSURANCE



Absolutely safe heating up to 300 °C thanks to the multi-level safety concept.



Intelligent control technology analyzes and reacts as required.



Intuitive operation:
Touch, turn & go!



Available as a Class 1 medical device (models UNm / UFm).

Model sizes: 30, 55, 75, 110, 160, 260, 450, 750, 1060

Model variants: SingleDISPLAY or TwinDISPLAY with natural convection (UN / UNm) or forced air circulation (UF / UFm).

Temperature setting ranges: +20 °C to +300 °C

All-round surface heating for optimal temperature distribution

The all-round surface heating is typical for Memmert and unique in the market. The protected heating elements of the four-sided all-round heating guarantee excellent temperature homogeneity.


Clever: Using the HeatBALANCE function enables the heating output, located in the lower and upper areas of the work area, to be regulated separately and be adaptable to any situation. This is done by the control of separate heating systems - directly on the device by the user.

Quality with added value for every facility

Memmert has always used stainless steel of the highest quality, because it is the ideal material for a heating oven. It is particularly corrosion-resistant, hygienic, robust and recyclable. Wherever Memmert devices are used, the superiority of high-quality stainless steel over painted sheet steel is evident, which can corrode very quickly after prolonged use.

Intuitive operation: Touch, turn & go!

The elegant ControlCOCKPIT, the operating unit with touch function, enables intuitive and uncomplicated control of parameters such as temperature, speed of the circulating air motor and the air flap position.

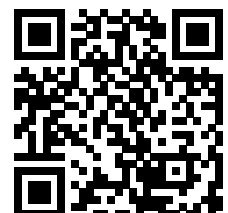



LABORATORY EQUIPMENT MAINTENANCE, REPAIR, CALIBRATION AND SALES

Established in 1987, Henderson Biomedical is the UK's leading laboratory equipment sales and service provider. Our knowledgeable team can provide you with excellent sales advice on a range of different types of laboratory equipment including centrifuges, refrigerators, freezers and heat sealers.

Henderson Biomedical is also able to provide you with first class after-sales service and calibration of your laboratory equipment. We are an **ISO 17025 (UKAS) accredited calibration test laboratory** and our team of Field Service Engineers cover the whole of the United Kingdom.

Please contact us for more information on the types of equipment we supply and the different after-sales services we can offer.

Henderson Biomedical

Unit 3, Swan Close
Croydon CR0 2DZ
United Kingdom

Tel: 020 8663 4610

**For sales enquiries: sales@henderson-biomedical.co.uk
For all other enquiries: info@henderson-biomedical.co.uk
www.henderson-biomedical.co.uk**

