


CTC256

Memmert environmental test chambers guarantee the perfect atmosphere for climate and temperature tests



On this page, you can find all the essential technical data on the Memmert environmental chamber CTC for climate testing and temperature testing. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at sales@memmert.com.

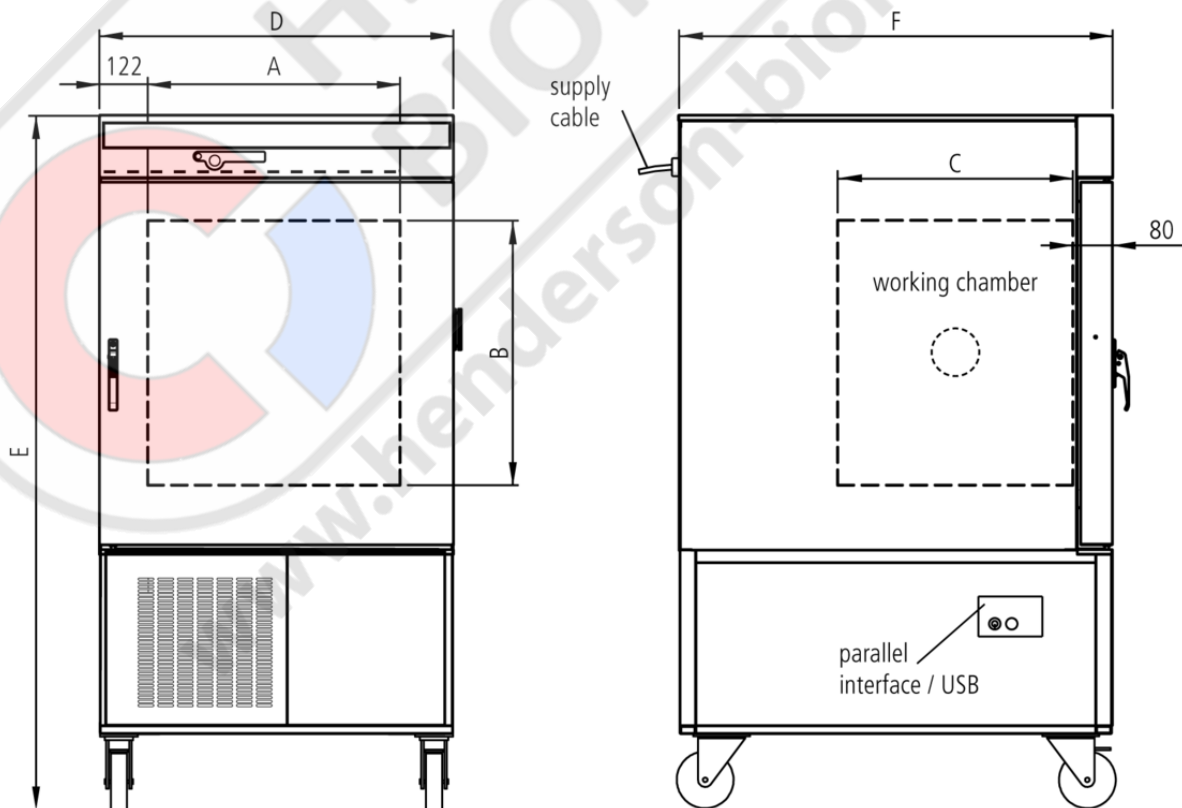

Temperature

Setting accuracy temperature	up to 99.9 °C: 0.1 / from 100 °C: 0.5
Temperature range	without humidity: from -42°C up to +190°C
Temperature range	with humidity from +10°C up to +95°C
Temperature distribution (spatial)	+/- 0.5 up to 2K
Temperature sensor	2 Pt100 sensors DIN Class A in 4-wire-circuit for mutual monitoring, taking over functions in case of an error
Display	resolution of display for setpoint values 0.1°C up to 99.9°C, 0.5°C from 100°C and for actual values 0.1°C (LED)

Humidity

Humidity control	active humidifying and de-humidifying adjustable from 10-98 % rh with digital display of relative humidity - resolution of display 0.5 %, setting accuracy 1 %
Humidification	

Control technology

Controller	Electronic microprocessor temperature controller with auto-diagnostic system
Set-up Display	digital display of all set parameters, such as temperature, weekdays, time, humidity, fan speed, programme status and set-up values - language to be chosen via set-up
Timer	integrated timer for tempering profiles of up to 40 ramps each, each segment adjustable from 1 min. to 999 hrs.
Calibration	three freely selectable temperature values, humidity 2-point calibration at 20% and 90 % rh

Ventilation

Fan	high-performance air fan in working chamber
Fan control	depending on operation status automatic adaption of fan speed resp. manual adjustment from 10 to 100 %

Communication

Interface USB	USB-interface incl. Memmert software "Celsius" for programming and documentation
Interface Printer	parallel printer interface (incl. real time clock with date function) for all PCL3-compatible ink jet printers for GLP-conforming documentation
Documentation	integrated ring memory as data logger for GLP-conforming long-term documentation of all relevant parameters - 1024 kB
Documentation	programme stored in case of power failure
Programming	chip-card control incl. 1 MEMORYCard XL with 32 kB storage capacity (max. 40 ramps)

Safety

Temperature control	mechanical temperature limiter TB, protection class 1 according to DIN 12880 to switch off the heating approx. 10°C above nominal temperature
Temperature control	over- and undertemperature monitor TWW, protection class 3.3
AutoSAFETY	additionally integrated over- and undertemperature protection "ASF", automatically following the setpoint value at a preset tolerance range, alarm in case of over- or undertemperature, heating is switched off in case of overtemperature, compressor in case of undertemperature
Autodiagnostic system	integral fault diagnostics for temperature and humidity control
Alarm	with visual and acoustic alarm in case of over-/under temperature and underhumidity, open door and empty water tank

Heating concept

Door	door heating to avoid condensation
Ring heater	high-performance ring heaters with optimised air circulation
Cooling	twin compressor; refrigerant R449A

Standard equipment

Works calibration certificate	standard values -20 °C and +160 °C, +30 °C at 60 % rh
Door	fully insulated stainless steel door with double-locking and 4-point adjustment, heated
Entry port	Entry port right, 80 mm, with stopper
Internals	1 stainless steel grid(s), electropolished
Interior	Easy-to-clean stainless steel interior, material 1.4301 (ASTM 304), hermetically welded

Stainless steel interior

Volume	256 l
Dimensions	$w_{(A)} \times h_{(B)} \times d_{(C)}$: 640 x 670 x 597 mm
Max. number of internals	6
Max. loading per internal	25 kg

Textured stainless steel casing

Dimensions	$w_{(D)} \times h_{(E)} \times d_{(F)}$: 898 x 1730 x 1100 mm (d +50mm door handle)
Installation	on lockable castors
Housing	rear zinc-plated steel

Electrical data

Electrical load 400 V, 3ph. w/o N 16A, 50 Hz / approx. 7000 W

Ambient conditions

Set Up The distance between the test chamber and the walls/ceiling must be at least 25 cm.

Ambient temperature 16 °C to 28 °C

Humidity rh max. 70 %, non-condensing

Altitude of installation max. 2,000 m above sea level

Pollution degree 2

Packing/shipping data

Transport information The appliances must be transported upright

Customs tariff number 8419 8998

Country of origin Federal Republic of Germany

WEEE-Reg.-No. DE 66812464


Dimensions approx incl. carton w x h x d: 1020 x 1910 x 1310 mm

Net weight approx. 337 kg

Gross weight carton approx. 463 kg

Temperature-humidity working range CTC

Not all climate chambers are the same. The humidity content of the chamber load, the ambient conditions and the respective temperature-humidity working range are decisive factors in the selection of the right appliance. In the adjacent diagram, you can see the possible temperature/humidity combinations for our climatic test chamber CTC. Within the respective temperature-humidity range, condensation-free permanent operation is possible. To which extent condensation may occur in the threshold range depends on the humidity content of the chamber load and the ambient conditions.

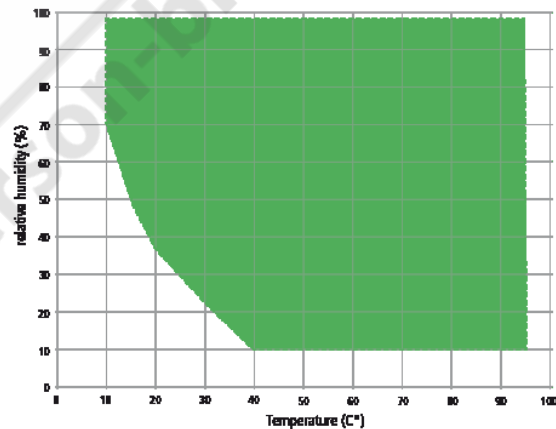

Temperature-humidity working range CTC

Standard units are safety-approved and bear the test marks

CE EAC





LABORATORY EQUIPMENT MAINTENANCE, REPAIR, CALIBRATION AND SALES

Established in 1987, Henderson Biomedical is the UK's leading laboratory equipment sales and service provider. Our knowledgeable team can provide you with excellent sales advice on a range of different types of laboratory equipment including centrifuges, refrigerators, freezers and heat sealers.

Henderson Biomedical is also able to provide you with first class after-sales service and calibration of your laboratory equipment. We are an **ISO 17025 (UKAS) accredited calibration test laboratory** and our team of Field Service Engineers cover the whole of the United Kingdom.

Please contact us for more information on the types of equipment we supply and the different after-sales services we can offer.

Henderson Biomedical

Unit 3, Swan Close
Croydon CR0 2DZ
United Kingdom

Tel: 020 8663 4610

For sales enquiries: sales@henderson-biomedical.co.uk
For all other enquiries: info@henderson-biomedical.co.uk
www.henderson-biomedical.co.uk

