

BioPlus EF600/660W



Biostorage you can depend on



BioPlus EF660W

The BioPlus EF600/660W is available in either white or stainless steel finish. All BioPlus models come with self-closing door.

Version 04 Valdid from 06.11.20

The **BioPlus EF600/660W** cabinet is a low temperature freezer that operates in the temperature range $-5/-35^{\circ}\text{C}$. It is a high capacity ventilated refrigeration unit that makes sure the interior quickly returns to the specified temperature after the door has been opened.

The control unit is a special Gram Multi Purpose Controller with LED display. It features acoustic alarm that provides clear notification if either maximum or minimum permitted temperature is reached. An acoustic door alarm is also fitted. All alarm parameters are set separately. Another embedded feature in the control-unit is voltage-free contact for IT networks, monitoring systems and/or remote alarms. Add to this alarm recording and Hi/Low temperature read-out. E-sensor for temperature reference in the storage compartment is standard.

The **BioPlus EF600/660W** utilises environmentally responsible HFC-free foam propellants and is available with HFC-free refrigerants.

The functional design ensures easy, ergonomically correct access to the storage space.



II 3G Ex nA nC nL IIB T2 Gc
TÜV 08 ATEX 354663 X

Biostorage
you can
depend on



BioPlus EF600/660W



Biostorage you can depend on

Top panel with BioLine control unit

The MPC unit is specially developed for users of BioLine storage cabinets. Among its many features are acoustic temperature and door alarms, and a voltage-free contact (NO/NC) for remote alarms. The display can be easily locked.



Access port

24.5 mm port for ease of access, e.g. sensors for external temperature surveillance.




E-sensor

Extra sensor for providing a temperature reference within the storage compartment.



Shelves and drawers

A wide range of versatile interior fittings and layouts is available. The operating specifications and the interior layout can both be customised.

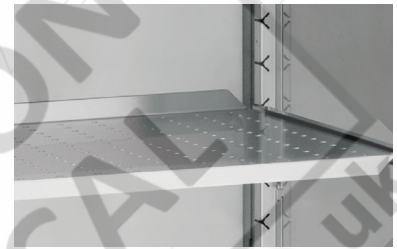


Inner doors

Draught shield that minimizes cold loss when the outer door is opened.

ATEX

All BioMidi models comply with ATEX Cat.3 Zone 2 requirements, both internal and externally, allowing for placement in locations where explosive gas atmospheres may be present.

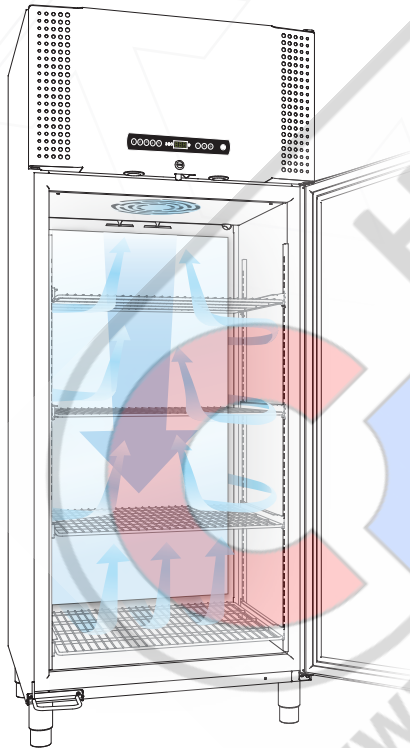

AIR DISTRIBUTION

The unique Gram BioLine air distribution system makes sure the temperature inside the cabinet remains stable at all times, while using the minimum amount of energy.

Technical specifications	EF600/660W
Temperature range	-5/-35 °C
Ambient temperature range	+10/+30 °C
Control Unit	Gram Control Unit with voltage-free contact, E-sensor, acoustic and visual door and temperature alarms that can be programmed individually, alarm recording and offset function
Material interior	Stainless steel and inner doors in PMMA
Material exterior	White lacquered steel or stainless steel
Dimensions mm (W x D x H)	600W: 837 x 756 x 1875/2125 / 660W: 837 x 756 x 2025/2275
Gross volume	600W: 600 litres/21.2 cubic feet / 660W: 660 litres/23.3 cubic feet
Net volume	600W: 432 litres/15.2 cubic feet / 660W: 484 litres/17.1 cubic feet
Net weight	Unpacked net weight without optional fittings: 600W: 146 kg / 660W: 152 kg
Modules for shelves and drawers	600W: Shelves 24 - Drawers 12 / 660W: Shelves 27 - Drawers 13
Insulation	60 mm polyurethane with HFC-free cyclopentane propellant
K-value	Solid door 0.3122/0.3131 W/m²K
Refrigerant	R290 / R404A/prepared for external compressor
CO ₂ e	R404A:902 / R404A (dual comp.): 1545
Refrigeration capacity at -40 °C	R290 244 watt; R404A 313 watt
Energy consumption	600W: 8.86 kWh/24h / 660W: 9.02 kWh/24h
Sound level	55,3 / 55,3 dB(A)
Base	Legs**
Connection	230 V, 50 Hz
Air system	Gram BioLine ventilated air distribution system
Defrost system	Automatic smart defrost with re-evaporation of defrost water

* Available in height less than 2000 mm

** The cabinets can also be delivered with castors or for mounting on a plinth at no extra charge





LABORATORY EQUIPMENT MAINTENANCE, REPAIR, CALIBRATION AND SALES

Established in 1987, Henderson Biomedical is the UK's leading laboratory equipment sales and service provider. Our knowledgeable team can provide you with excellent sales advice on a range of different types of laboratory equipment including centrifuges, refrigerators, freezers and heat sealers.

Henderson Biomedical is also able to provide you with first class after-sales service and calibration of your laboratory equipment. We are an **ISO 17025 (UKAS) accredited calibration test laboratory** and our team of Field Service Engineers cover the whole of the United Kingdom.

Please contact us for more information on the types of equipment we supply and the different after-sales services we can offer.

Henderson Biomedical

Unit 3, Swan Close
Croydon CR0 2DZ
United Kingdom

Tel: 020 8663 4610

**For sales enquiries: sales@henderson-biomedical.co.uk
For all other enquiries: info@henderson-biomedical.co.uk
www.henderson-biomedical.co.uk**

