

# BioPlus 500



Biostorage you can depend on




The **BioPlus 500** range is designed for the storage of the most delicate biomaterials, in situations where even tiny fluctuations in conditions inside the storage cabinet can have a serious effect on the contents.

BioPlus cabinets also enable you to reduce internal relative humidity, thus cutting down on any likelihood of undesirable contaminants getting near delicate biomaterials.

In terms of performance, these units are designed to provide the very best results even under exceptional conditions. The BioPlus utilises environmentally responsible HFC-free foam propellants and is available with HFC-free refrigerants.

BioPlus 500.  
Drawers and glass door are optional extras.



 II 3G Ex nA nC nL IIB T2 Gc  
TÜV 08 ATEX 354663 X

The BioPlus 500 is available as a refrigerator or freezer in either white or stainless steel finish. All BioPlus models come with self-closing door. Glass door is available for refrigerator.

Biostorage  
you can  
depend on

Version 03. Valid from 15.10.20



# BioPlus 500



Biostorage you can depend on

## Top panel with BioLine control unit

The MPC unit is specially developed for users of BioLine storage cabinets. Among its many features are acoustic temperature and door alarms, and a voltage-free contact (NO/NC) for remote alarms. The display can be easily locked.



## Access port

24.5 mm port for ease of access, e.g. sensors for external temperature surveillance.



## E-sensor

Extra sensor for providing a temperature reference within the storage compartment.



## Shelves and drawers

A wide range of versatile interior fittings and layouts is available. The operating specifications and the interior layout can both be customised.




## Stainless steel interior

The internal sides and bottom as well as shelf and wall rails and support rails are made of stainless steel.

## ATEX

All BioPlus models comply with ATEX Cat.3 Zone 2 requirements, both internal and externally, allowing for placement in locations where explosive gas atmospheres may be present.

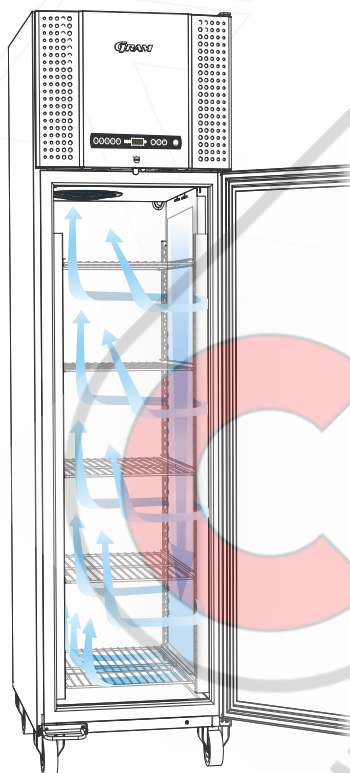

## AIR DISTRIBUTION

The unique Gram BioLine air distribution system makes sure the temperature inside the cabinet remains stable at all times, while using the minimum amount of energy.

Technical Specifications	BioPlus ER500	BioPlus RF500
Temperature range	+20/-2 °C	-5/-25 °C
Ambient temperature range	Solid door +10/+43 °C; glass door (ER) +10/+38 °C	
Control Unit	Gram Control Unit with voltage-free contact, E-sensor, acoustic and visual door and temperature alarms that can be programmed individually, alarm recording and offset function.	
Material interior	Stainless steel	
Material exterior	White lacquered steel or stainless steel	
Dimensions mm (W x D x H)	600 x 806 x 2025/2275	
Gross volume	500 litres/17.7 cubic feet	
Net volume	365 litres/12.9 cubic feet	
Net weight	Unpacked net weight without optional fittings: 130 kg	
Modules for shelves and drawers	Shelves 52 - Drawers 13	
Insulation	60 mm polyurethane with HFC-free cyclopentane propellant	
K-Value	Glass door 0.4641 W/m²K	Solid door 0.3149 W/m²K
Refrigerant	R290 / R134a	R290 / R404A
Co <sub>2</sub> e	R134a: 372	R404A: 981
Refrigeration capacity at -10 °C	R290 389 Watt / R134a 314 Watt	-
Refrigeration capacity at -25 °C	-	R290 374 Watt / R404a 616 Watt
Energy consumption	1.26 kWh/24h	3.83 kWh/24h
Sound level	46,2 dB(A)	48,7 dB(A)
Base	Legs*	
Connection	230 V, 50 Hz	
Air system	Gram BioLine ventilated air distribution system	
Defrost system	Automatic smart defrost with re-evaporation of defrost water	

\* The cabinets can also be delivered with castors or for mounting on a plinth at no extra charge.



***LABORATORY EQUIPMENT MAINTENANCE, REPAIR, CALIBRATION AND SALES***

Established in 1987, Henderson Biomedical is the UK's leading laboratory equipment sales and service provider. Our knowledgeable team can provide you with excellent sales advice on a range of different types of laboratory equipment including centrifuges, refrigerators, freezers and heat sealers.

Henderson Biomedical is also able to provide you with first class after-sales service and calibration of your laboratory equipment. We are an **ISO 17025 (UKAS) accredited calibration test laboratory** and our team of Field Service Engineers cover the whole of the United Kingdom.

Please contact us for more information on the types of equipment we supply and the different after-sales services we can offer.

**Henderson Biomedical**

Unit 3, Swan Close  
Croydon CR0 2DZ  
United Kingdom

**Tel: 020 8663 4610**

**For sales enquiries: [sales@henderson-biomedical.co.uk](mailto:sales@henderson-biomedical.co.uk)  
For all other enquiries: [info@henderson-biomedical.co.uk](mailto:info@henderson-biomedical.co.uk)  
[www.henderson-biomedical.co.uk](http://www.henderson-biomedical.co.uk)**

