

BioPlus 930



Biostorage you can depend on



BioPlus 930



The BioPlus 930 is available as a refrigerator or freezer in either white or stainless steel finish. All BioPlus models come with self-closing door. Glass door is available for refrigerator.

The **BioPlus 930** range is designed for the storage of the most delicate biomaterials, in situations where even tiny fluctuations in conditions inside the storage cabinet can have a serious effect on the contents.

BioPlus ER cabinets also enable you to reduce internal relative humidity, thus cutting down on any likelihood of undesirable contaminants getting near delicate biomaterials.

In terms of performance, these units are designed to provide the very best results even under exceptional conditions. The BioPlus utilises environmentally responsible HFC-free foam propellants and is available with HFC-free refrigerants.



II 3G Ex nA nC nL IIB T2 Gc

Biostorage
you can
depend on

Version 03. Valid from 19.10.20



BioPlus 930



Biostorage you can depend on

Top panel with BioLine control unit

The MPC unit is specially developed for users of BioLine storage cabinets. Among its many features are acoustic temperature and door alarms, and a voltage-free contact (NO/NC) for remote alarms. The display can be easily locked.



Access port

24.5 mm port for ease of access, e.g. sensors for external temperature surveillance.



E-sensor

Extra sensor for providing a temperature reference within the storage compartment.



Shelves and drawers

A wide range of versatile interior fittings and layouts is available. The operating specifications and the interior layout can both be customised.




Stainless steel interior

The internal sides and bottom as well as shelf and drawer brackets and support rails are made of stainless steel.

ATEX

All BioPlus models comply with ATEX Cat.3 Zone 2 requirements, both internal and externally, allowing for placement in locations where explosive gas atmospheres may be present.

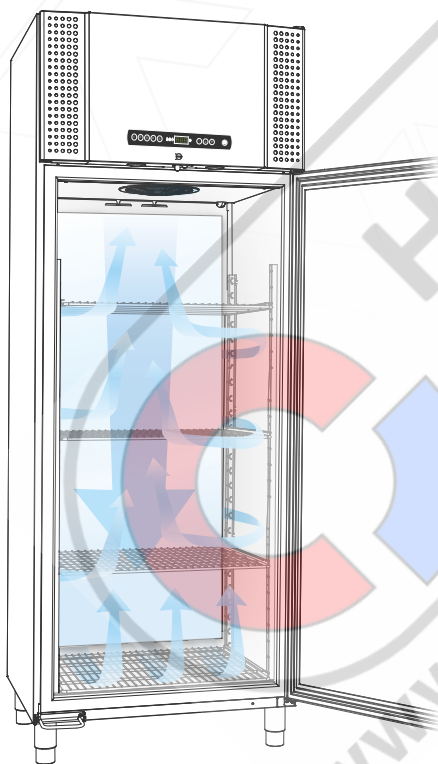

AIR DISTRIBUTION

The unique Gram BioLine air distribution system makes sure the temperature inside the cabinet remains stable at all times, while using the minimum amount of energy.

Technical specifications	BioPlus ER930	BioPlus RF930
Temperature range	+20/-2 °C	-5/-25 °C
Ambient temperature range	Solid door +10/+43 °C; Glass door +10/+38 °C	
Control Unit	Gram Control Unit with voltage-free contact, E-sensor, acoustic and visual door and temperature alarms that can be programmed individually, alarm recording and offset function.	
Material interior	Stainless steel	
Material exterior	White lacquered steel or stainless steel	
Dimensions mm (W x D x H)	780 x 1045 x 2025/2275	
Gross volume	930 litres/32.8 cubic feet	
Net volume	702 litres/30.0 cubic feet	
Net weight	Unpacked net weight without optional fittings: 165 kg	
Modules for shelves and drawers	Shelves 27	
Insulation	60 mm polyurethane with HFC-free cyclopentane propellant	
K-Value	Glass door 0.4563 W/m²K	Solid door 0.3137 W/m²K
Refrigerant	R290 / R134a	R290 / R404A
Co ₂ e	R134a: 572	R404A: 1177
Refrigeration capacity at -10 °C	R290 389 watt; R134a 314 watt	-
Refrigeration capacity at -25 °C	-	R290 627 watt; R404A 487 watt
Base	Legs*	
Energy consumption	2.10 kWh/24h	6,12 kWh/24h
Sound level	48,5 db(A)	49,5 db(A)
Connection	230 V, 50 Hz	
Air system	Gram BioLine ventilated air distribution system	
Defrost system	Automatic smart defrost with re-evaporation of defrost water	

* The cabinets can also be delivered with castors or for mounting on a plinth at no extra charge



LABORATORY EQUIPMENT MAINTENANCE, REPAIR, CALIBRATION AND SALES

Established in 1987, Henderson Biomedical is the UK's leading laboratory equipment sales and service provider. Our knowledgeable team can provide you with excellent sales advice on a range of different types of laboratory equipment including centrifuges, refrigerators, freezers and heat sealers.

Henderson Biomedical is also able to provide you with first class after-sales service and calibration of your laboratory equipment. We are an **ISO 17025 (UKAS) accredited calibration test laboratory** and our team of Field Service Engineers cover the whole of the United Kingdom.

Please contact us for more information on the types of equipment we supply and the different after-sales services we can offer.

Henderson Biomedical

Unit 3, Swan Close
Croydon CR0 2DZ
United Kingdom

Tel: 020 8663 4610

**For sales enquiries: sales@henderson-biomedical.co.uk
For all other enquiries: info@henderson-biomedical.co.uk
www.henderson-biomedical.co.uk**

