




RAFR21043 Solid Door Advanced Laboratory Refrigerator

Part Number	RAFR21043
Capacity (Litres)	650
Style	Freestanding
Exterior Dimensions (HxWxD)	1985 x 710 x 800mm
Energy Consumption	2.43 kWh/24h
Operating Temperature	+5°C
Temperature Set Point Range	+2°C to +10°C
Interior Construction	Stainless Steel
Exterior Construction	White polyester coated sheet steel
Shelves	4
Maximum shelf loading kg	70
Shelf dimensions (W x D)	500 x 620 mm
Door Hinging	Right side hinging only
Door Style	Solid
Internal Light	✓
Lockable Castors	✓
Door Lock	✓
Access Port	Optional Extra
Digital Temperature Control	✓
Digital Temperature Display	✓
High Temperature Alarm	✓
Low Temperature Alarm	✓
Door Open Alarm	✓
Maximum Temperature Display	✓
Minimum Temperature Display	✓
Auto Defrost	✓
Fan Assisted Cooling	✓
Low GWP Refrigerant	✓
Ammonia Free	✓

Robustly constructed from steel with a stainless steel interior this refrigerator from the **Labcold** Advanced Laboratory Range is designed to provide close temperature control for temperature sensitive materials in the laboratory.

In addition, remote alarm contacts can be fitted as well as access ports. Contact us for bespoke racking inventory systems or specialised applications.

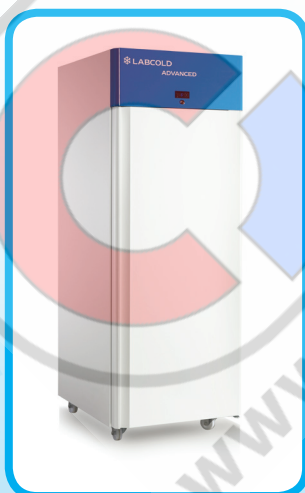

The air is continually circulated internally to ensure that where ever product is placed inside the cabinet, it stays within range. There is no internal 'step' and the interior is hygienic stainless steel making it easy to clean and stock. The 4 heavy duty shelves are fully adjustable and more shelves can be purchased as an optional extra. The robust doors are independently lockable and for increased security and convenience are self closing.

The cabinet is fitted with heavy duty, lockable castors making it easy to manoeuvre and clean behind which means this Advanced Laboratory Refrigerator is ideal for demanding environments.



LABORATORY EQUIPMENT MAINTENANCE, REPAIR, CALIBRATION AND SALES

Established in 1987, Henderson Biomedical is the UK's leading laboratory equipment sales and service provider. Our knowledgeable team can provide you with excellent sales advice on a range of different types of laboratory equipment including centrifuges, refrigerators, freezers and heat sealers.

Henderson Biomedical is also able to provide you with first class after-sales service and calibration of your laboratory equipment. We are an **ISO 17025 (UKAS) accredited calibration test laboratory** and our team of Field Service Engineers cover the whole of the United Kingdom.

Please contact us for more information on the types of equipment we supply and the different after-sales services we can offer.

Henderson Biomedical

Unit 3, Swan Close
Croydon CR0 2DZ
United Kingdom

Tel: 020 8663 4610

For sales enquiries: sales@henderson-biomedical.co.uk
For all other enquiries: info@henderson-biomedical.co.uk
www.henderson-biomedical.co.uk

